


RT-SPOT[®]


Beekley
Medical[®]
when your diagnosis
must be right[®]

The only skin marker for CT simulation that is elevated from the skin line and helps to exclude the marker from the automatic body contour


- ▶ Minimizes dose perturbation and helps ensure dose accuracy
- ▶ Saves time by reducing the need for manual contouring and density override
- ▶ Flexible and easier to form into desired shapes
- ▶ Stays in place on the skin and does not lift
- ▶ Lower density, non-metallic line, minimizes artifact

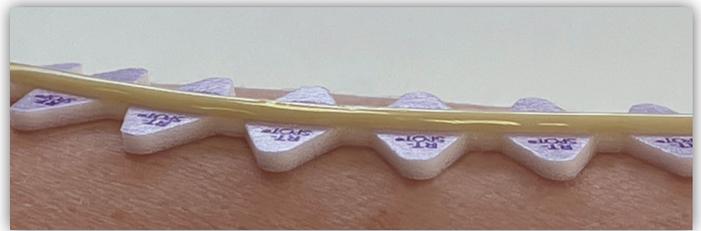

▶ Medical grade adhesive foam backing elevates RT-SPOT from the skin line

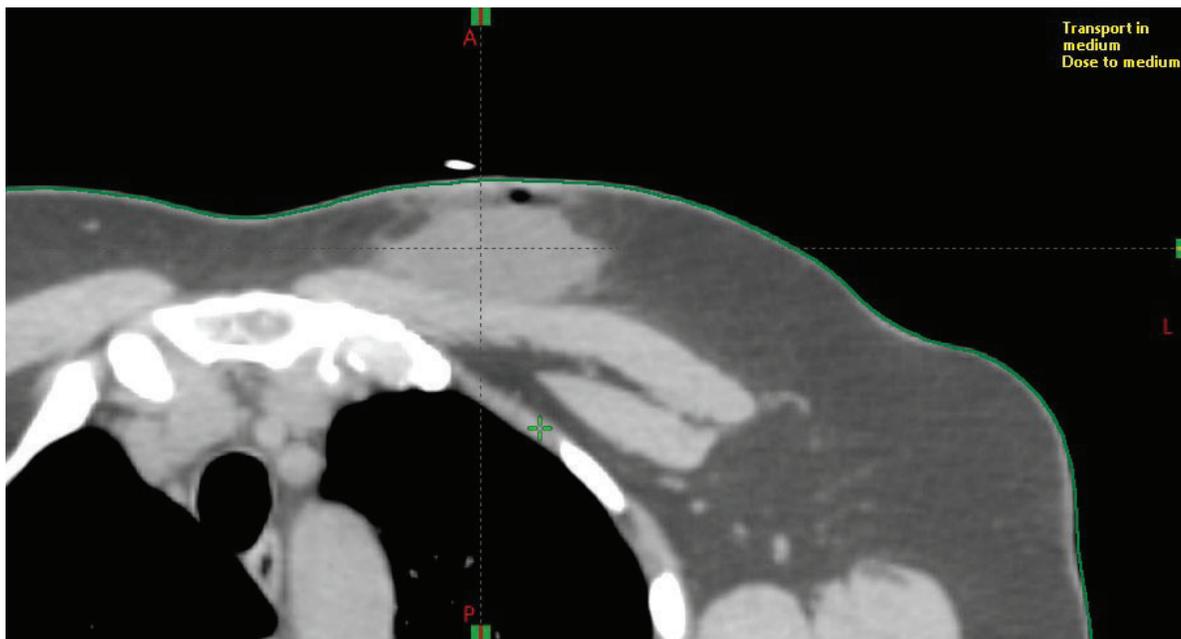

▶ Axial CT Scan of female breast radiation patient with surgical scar identified by RT-SPOT skin marker for CT simulation. RT-SPOT excluded from body contour.

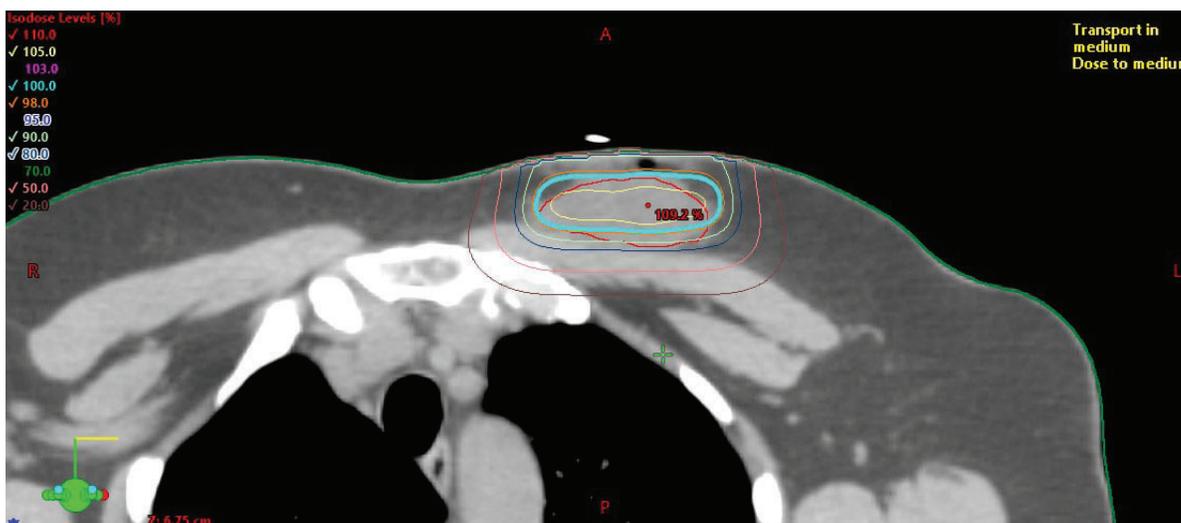

▶ Corresponding image of electron dose distribution unperturbed by RT-SPOT skin marker for CT simulation.

RT-SPOT[®]

Beekley
Medical[®]
when your diagnosis
must be right[®]

Existing skin markers affect dose calculations and must be manually contoured out, increasing the dose planning time per patient.

Manual contouring can be imperfect leading to dose inaccuracies or rework of treatment plans.

RT-SPOT is an elevated skin marker designed to minimize dose perturbation and reduce the need for manual contouring.

RT-SPOT skin marker for CT simulation is a flexible line that contours more easily into desired shapes and stays in place without lifting. Ideal for marking field borders, tangents, scars, sarcomas and larger treatment areas.



▲ The moldable backing allows RT-SPOT skin marker for CT simulation to bend into various shapes – even small circles



REF 127 QTY 140cm cut to measure roll

- 2.0mm diameter, elevated, non-metallic line
- For marking field borders, tangents, and scars

1.800.233.5539 • info@beekley.com • beekley.com

Customers outside the U.S. – Contact your local distributor for pricing and product availability.

To locate a distributor, call +1.860.583.4700 or email international@beekley.com

 Manufactured by Beekley Corporation
One Prestige Lane, Bristol, CT 06010-7454 USA
Tel: 1.800.233.5539 or +1.860.583.4700
Fax: 1.800.735.1234 beekley.com

BEEKLEY, BEEKLEY MEDICAL, RT-SPOT, SPOT, and WHEN YOUR DIAGNOSIS MUST BE RIGHT are Reg. U.S. Pat. & Tm. Off. BEEKLEY, RT-SPOT, and SPOT are Registered Community Trademarks. BEEKLEY is a registered trademark in Canada. © 2023 Beekley Corporation. All rights reserved.

WARNINGS AND PRECAUTIONS:

- Do Not Use in Skin Folds
- Do Not Use Under Physical Bolus
- Single Patient Use Only; Do Not Reuse 
- External Use Only; Do Not Ingest or Swallow
- Do Not Place on Compromised Skin
- Remove After Use
- Not Made With Natural Rubber Latex
- Lead-Free
- Non-Sterile 

   EMERGO EUROPE
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands



REV: RTS_INF_0523